




LIFECOMFORT® HYBRID COMFORT MATTRESS

The Hybrid Comfort is an Acute mattress system offering foam-in-cell technology enabling flexible operation as a high comfort foam static support surface or a powered-up alternating pressure support surface. Specifically re-engineered for increased patient comfort, foam layering now effectively supports patient weight ranging from 30 to 380kg.

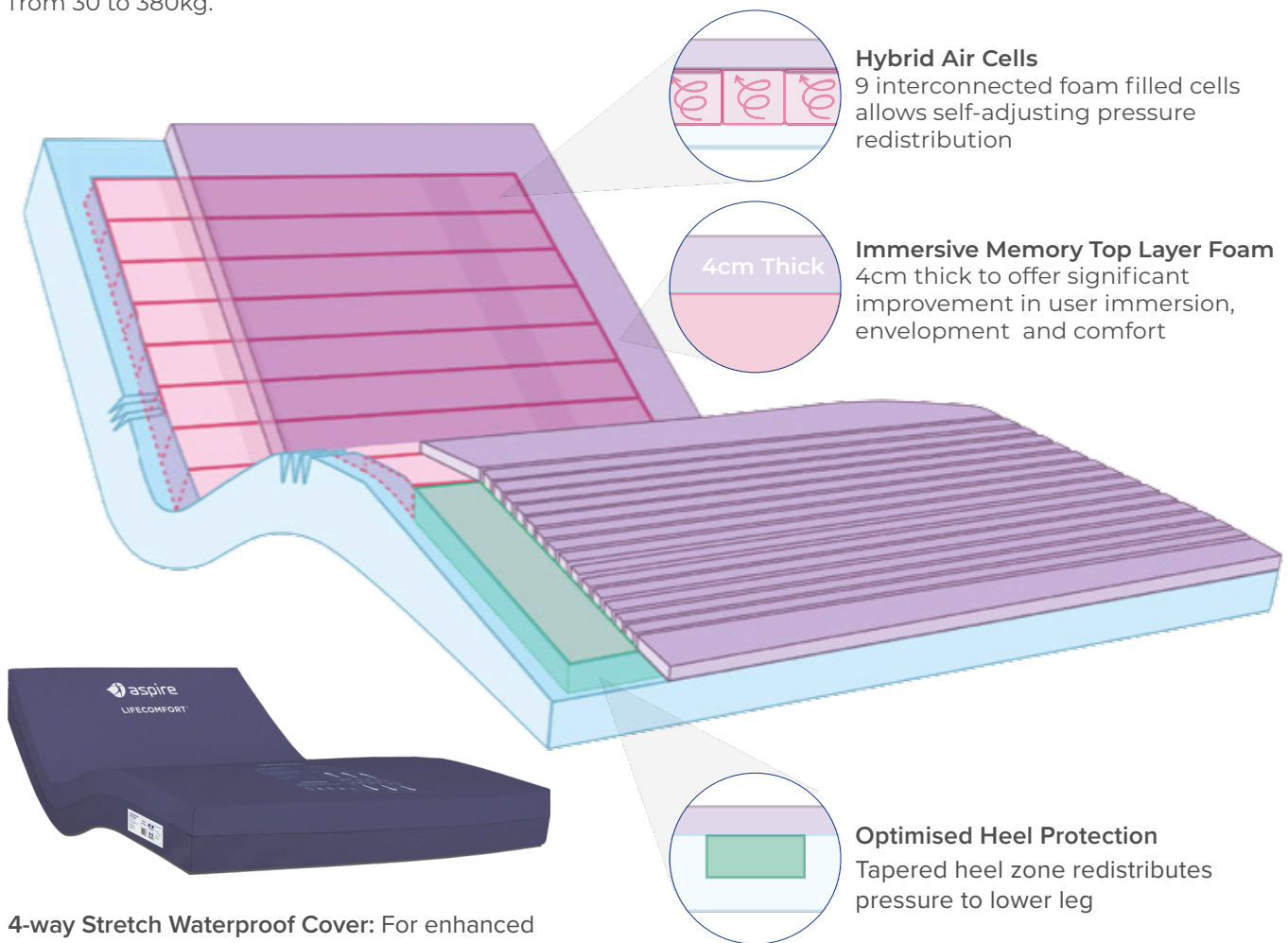

4-way Stretch Waterproof Cover: For enhanced shear mitigation and friction reduction

Intelligent Hose Management: When not in use the air-hose can be stored internally, reducing infection control and OH&S risks

Integrated Shear Reduction: Memory foam layer moves independently of the cell surface reducing shear when bed is adjusted

Innovative Pump Technology: Reduces noise and vibration for undisturbed sleep



Digital Pump sold separately
Single and King Single - PTM638030*
Bariatric - PTM638032**

APPLICATIONS

- Aged Care
- Hospital
- Ambulant patients
- Palliative Care
- Home Care




5 YEAR WARRANTY
Cover and Foam





Asset Management Enabled through unique batch serial numbering

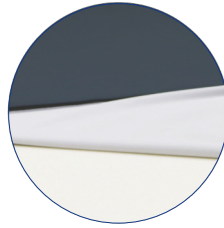

Hospital-grade Antimicrobial Foam Ensures superior infection prevention



Infection Control Fully welded seams with waterfall zip covers, protect against fluid ingress

CLINICAL BENEFITS

Reduced Pressure Injury Risk:

- Medium-High mattress risk rating with 4-way stretch cover providing high MVTR for optimal skin microclimate
- Memory foam layer with integrated shear reduction moves independently during bed adjustment
- Tapered heel zone redistributes pressure; castellated surface reduces shear through wave-like movement
- Enhanced visco-elastic memory foam surface (doubled thickness) improves immersion and envelopment
- Self-adjusting pressure redistribution via 9 interconnected foam cells allowing free air movement
- Hybrid system enables seamless conversion from static to alternating pressure support

Reduced Clinical Impact:

- Prevents multidisciplinary intervention requirements (medical, nursing, dietetics, AHP)
- Eliminates extended offloading periods (weeks to months for healing)
- Reduces increased support needs: personal care, domestic assistance, meal preparation

Enhanced Comfort:

- Thickened visco-elastic memory foam improves user immersion and comfort
- Integrated memory foam overlay with self-adjusting pressure redistribution
- Innovative pump technology minimizes noise/vibration for uninterrupted sleep

Improved Mobility:

- Firm foam bed edges enhance transfer safety and independence



Clinical Assistance Guide

LIFE CYCLE & COMPLIANCE

Warranty: 5 years Foam | 5 years Cover

Machine Washable Covers: Up to 95°

Durability + Performance: ISO 3385

Safe Working Load Verified: AS/NZ 2282

Flame Retardant Foam + Covers: AS/NZ 4088.1 Cribs 5

TGA: Lifecomfort Mattresses are registered in the ARTG, AUST R ARTG 508789

Code	Description	Length	Width	Height
PTM638005	Lifecomfort Hybrid Comfort - Single	1980	895	170
PTM638006	Lifecomfort Hybrid Comfort - King Single	1980	1050	170
PTM638007**	Lifecomfort Hybrid Comfort - Bariatric	1980	1050	170

DVA AB14 *AB19 NDIS 05_043306211_0103_1_2 *05_043306211_0103_1_2 **05_043306804_0103_1_2